

Angio Chest Flash for Dissection

Siemens Flash

Application Examples: aneurysm, stenosis, vascular anomalies

Oral Contrast	1 glass H2O
IV Contrast / Volume	Omnipaque 350 / MedRad P3T
Injection Rate	P3T / 3.5 – 7 mL/sec

Detector Collimator	Acq 128 x 0.6 mm
Care kV	On / 120 kV
Care Dose 4D	On / 65 mAs
Rotation Time (seconds)	0.28
Pitch	0.6
Typical CTDIvol	4.40 mGy \pm 50%

Technical Factors

Care Bolus ROI Location / HU	Aortic Arch / 180
Monitoring Delay	10 seconds
Cycle Time	1.43 seconds
Scan Delay	3 seconds
Breath Hold	No

Detector Collimator	Acq 128 x 0.6 mm
Care kV	On / 120 kV
Care Dose 4D	On / 130 mAs
Rotation Time (seconds)	0.28
Pitch	2.2
Typical CTDIvol	6.23 mGy \pm 50%

Topogram: Lateral and AP, 512 mm

Chest w/o	Recon Type	Width/Increment	Algorithm	Safire	Window	Series Description	Networking	Post Processing
Recon 1	Axial	5 x 5	I41f	2	Mediastinum	AXIAL WITHOUT	PACS	None


AngioChest	Recon Type	Width/Increment	Algorithm	Safire	Window	Series Description	Networking	Post Processing
Recon 1	Axial	3 x 3	I40f	2	Mediastinum	AXIAL	PACS	None
Recon 2	3D:COR	3 x 3	I30f	2	Angio	COR MIP	PACS	Coronal MIP
Recon 3	3D:SAG	3 x 3	I30f	2	Angio	SAG MIP	PACS	Sagittal MIP
Recon 4	3D:OBL	3 x 3	I30f	2	Angio	OBL MIP	PACS	Oblique MIP
Recon 5	Axial	0.6 x 0.6	I26f	2	Angio	AXIAL 0.6 STND	TR & PACS	None

IV Placement: ≥ 18 gauge, *preferably* in antecubital (AC) fossa.**Patient Preparation:** Give one glass of water just prior to scan.**Patient Position:** Patient lying supine with arms comfortably above head and lower legs supported.**Scan Instructions:** Place pre-monitoring ROI in aortic arch, away from any calcium deposits.**Scan Range:** Lung apices through celiac artery.**Recons and Reformats:** Adjust FoV to fit body contour.


3x3 Coronal MIP



3x3 Sagittal MIP

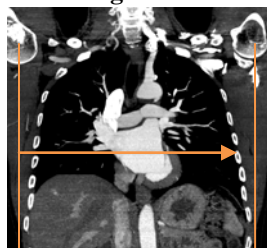

3x3 Oblique MIP

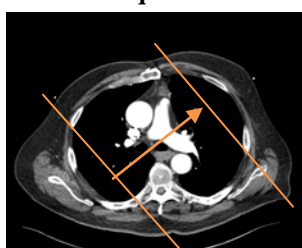**3D:** Upon request.