


Post Treatment Head / Brain Mets - Mulithance

Application Examples: post tumor removal, post radiation therapy


Sequence	FoV	Thickness	Gap
SAGITTAL T1	24CM	4MM	1MM
AX T2 TSE PROPELLER	24CM	4MM	1MM
AX T2 FLAIR PROPELLER	24CM	4MM	1MM
AX SWAN	24CM	4MM	1MM
AX T1	24CM	4MM	1MM
AX DWI (B-0, B-1000)	24CM	3MM	1MM
INJECT CONTRAST (1mmol/Kg)			
SAG CUBE FLAIR (<i>RECON AX & COR 1 x1</i>)	24CM	1.2MM	NO GAP
AX T1 FS SPGR WITH GAD	24CM	4MM	1MM
AX BRAVO WITH GAD (<i>Recon in Sagittal and Cor 1x1</i>)	25.6CM	2MM	NO GAP



AXIAL Coverage

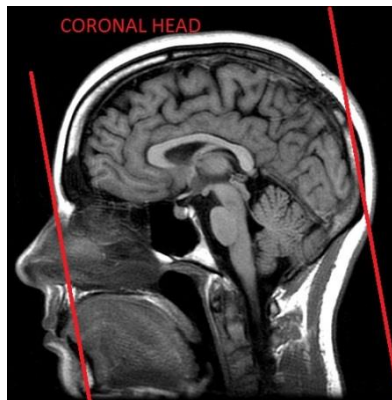

COR Coverage

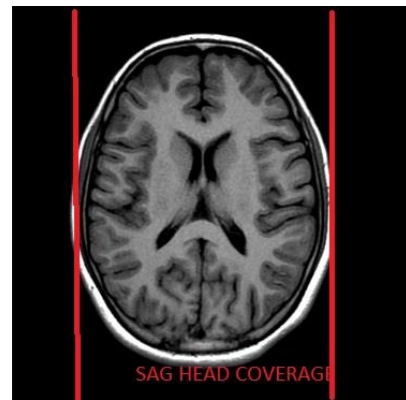

SAG Coverage