


## ***IAC / POST OP CHOLESTEATOMA***


Application Examples: hearing loss, follow up schwannoma, vertigo

Coil to Use: Head Matrix

Sequence	FoV	Thickness/Spacing	Gap	Nex/Aver	Matrix
AX T2 FLAIR PROPELLER ( <i>Whole Brain</i> )	24CM	4MM	1MM	1.5	288/288
AX DWI PROPELLER (B-0, B-1000) <b>Send Prop Combo to PACS</b>	24CM	3MM	NO GAP	2	128/128
AX T1 FSE IAC'S	18CM	3MM	NO GAP	3	288/192
COR T1 FSE IAC'S	18CM	3MM	NO GAP	3	288/192
AX FIESTA IAC'S ( <b>RECON COR 1X1</b> )	18CM	1MM	NO GAP	1	256/320
INJECT CONTRAST (1mmol/Kg)					
COR T1 FS SPGR IAC'S WITH GAD	18CM	3MM	NO GAP	2	320/224
AXIAL T1FS SPGR IAC'S WITH GAD	18CM	3MM	NO GAP	2	320/224
SAG FSPGR BRAVO WITH GAD <b>(RECON AX AND COR 1x1) whole heads</b>	24CM	2MM	NO GAP	1	256/256

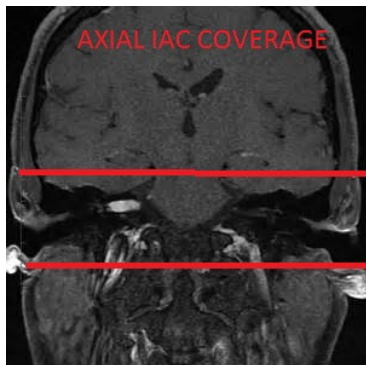

AXIAL Coverage

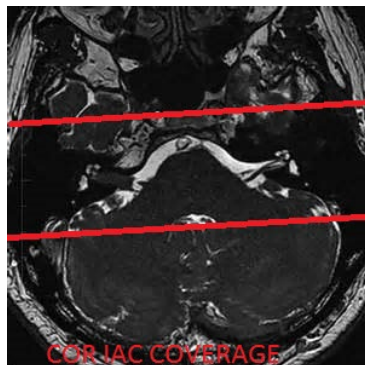

COR Coverage