



Pituitary

Application Examples: follow up micro adenoma or macro adenoma, growth hormone deficiency, low testosterone levels, acromegaly

Sequence	FoV	Thickness	Gap
COR T1	15CM	3MM	NO GAP
COR T2	15CM	3MM	NO GAP
SAG T1	15CM	3MM	NO GAP
SAG T2	15CM	3MM	NO GAP
INJECT CONTRAST (1mmol/Kg)			
COR DYNAMIC WITH GAD	18CM	2.5MM	NO GAP
COR T1 FS WITH GAD	15CM	3MM	NO GAP
SAG T1 FS WITH GAD	24CM	2MM	NO GAP

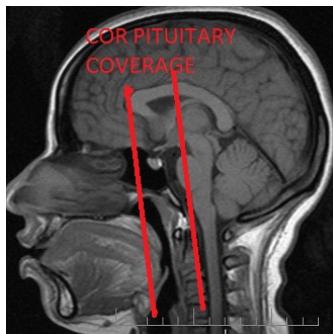

COR Coverage

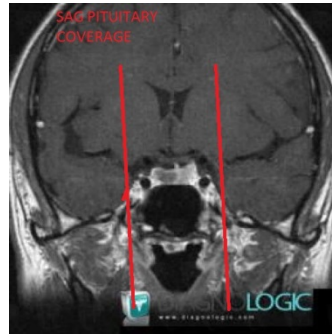

SAG Coverage