


Orbits

Application Examples: vision loss, optic neuritis

Coverage: Axial coverage is top orbits to mid maxillary sinus and coronal coverage is anterior eye thru sella.

Sequence	FoV	Thickness	Gap
COR T1 ORBITS	18CM	3MM	NO GAP
COR T2 FS ORBITS	18CM	3MM	NO GAP
COR T2 FLAIR FS ORBITS	18CM	3MM	NO GAP
AX FIESTA ORBITS <i>(RECON COR 1x1)</i>	18CM	0.6MM	NO GAP
AX T2 FS ORBITS	18CM	3MM	NO GAP
AX DWI WHOLE HEAD <i>(VALUES B-0,B-1000)</i>	24CM	3MM	1MM
INJECT CONTRAST (1mmol/Kg)			
COR T1 FS SPGR ORBITS WITH GAD	18CM	3MM	NO GAP
AX T1 FS SPGR ORBITS WITH GAD	18CM	3MM	NO GAP

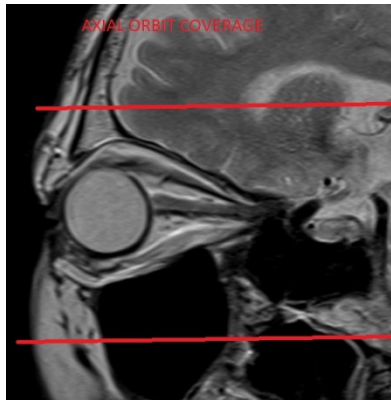

AXIAL Coverage



COR Coverage