


MS Head

Application Examples: facial weakness, history of multiple sclerosis, arm or leg weakness in patient under 40

Sequence	FoV	Thickness	Gap
SAGITTAL SPGR NEUROQUANT- RECON AX/SAG ONLY 1x1 Send neuroquant sequence to lesionquant	25.6CM	1MM	NO GAP
AX T2 PROPELLER	24CM	4MM	1MM
AX T2 FLAIR LESIONQUANT – SEND TO LESIONQUANT	23CM	3MM	1MM
AX SWAN	24CM	4MM	1MM
AX DWI (B-0, B-1000)	24CM	3MM	1MM
INJECT CONTRAST (1mmol/Kg)			
SAG CUBE FLAIR WITH GAD (RECON AX & COR 1 x1)	24CM	1.2MM	NO GAP
SAG T1 FS SPGR WITH GAD	24CM	4MM	1MM
AX T1 FS SPGR WITH GAD	24CM	4MM	1MM

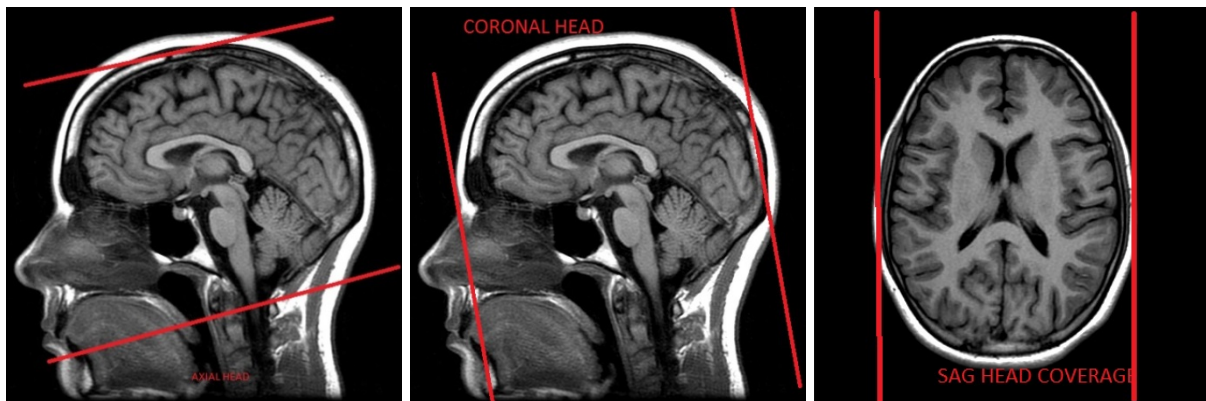

AXIAL Coverage

COR Coverage

SAG Coverage