


Skull Base


Application Examples: trigeminal neuralgia, oral cancer

Coverage: axial coverage is top orbits to C4 and coronal coverage is anterior orbits thru back of head.

Sequence	FoV	Thickness	Gap
COR T1 SKULL BASE	18CM	3MM	NO GAP
COR T2 FS SKULL BASE	18CM	3MM	NO GAP
AX T2 FS SKULL BASE	18CM	3MM	NO GAP
AX T1 SKULL BASE	18CM	3MM	NO GAP
AX DWI WHOLE HEAD TO INCLUDE ENTIRE SKULL BASE (VALUES B-0, B-1000)	24CM	3MM	1MM
INJECT CONTRAST (1mmol/kg)			
COR T1 FS SPGR SKULL BASE WITH GAD	18CM	3MM	NO GAP
AX T1 FS SPGR SKULL BASE WITH GAD	18CM	3MM	NO GAP

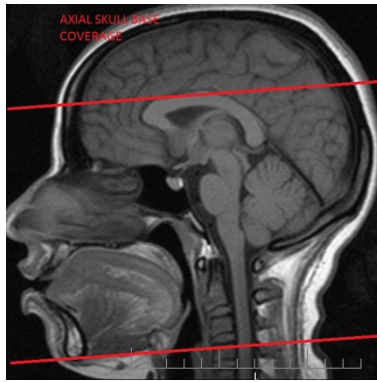

AXIAL Coverage

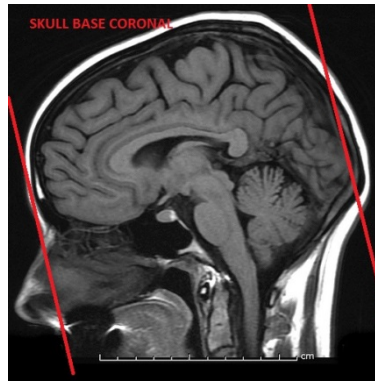

COR Coverage