



Cervical Spine

Application Examples: DJD, neck pain, arm pain


Sequence	FoV	Thickness	Gap
SAG T1	24CM	3MM	0.5MM
SAG T2	24CM	3MM	0.5MM
SAG STIR	24CM	3MM	0.5MM
AX 2D MERGE / MEDIC (SIEMENS)	18CM	3MM	0.3MM
AX T2	18CM	3MM	0.3MM
RIGHT SAG OBL T2	24CM	2MM	0.5MM
LEFT SAG OBL T2	24CM	2MM	0.5MM

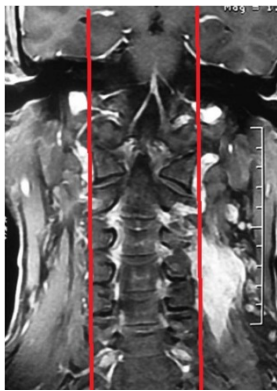

SAG Coverage

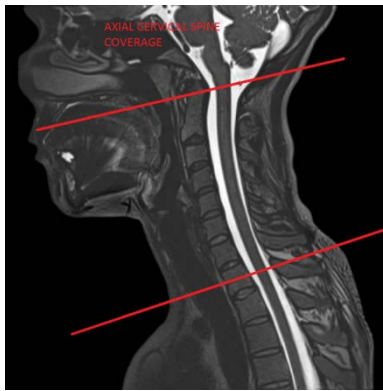

AXIAL Coverage

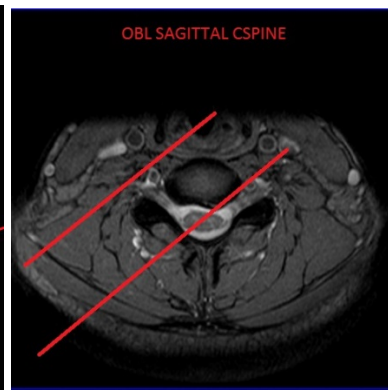

OBL T2 Coverage