



Cervical Spine MS / Tumor / Infection

Application Examples: MS, tumor, infection

Sequence	FoV	Thickness	Gap
SAG T1	24CM	3MM	0.5MM
SAG T2	24CM	3MM	0.5MM
SAG STIR	24CM	3MM	0.5MM
AX 2D MERGE / MEDIC (SIEMENS)	18CM	3MM	0.3MM
AX T2	18CM	3MM	0.3MM
AX DWI (<i>VALUES B-0, B-1000</i>)	24CM	3MM	0.5MM
INJECT CONTRAST (1mmol/Kg)			
SAG T1 FS WITH GAD	24CM	3MM	0.5MM
AX T1 FS WITH GAD	18CM	3MM	0.3MM

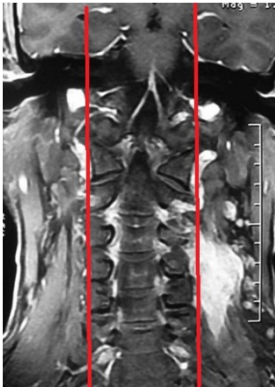

SAG Coverage

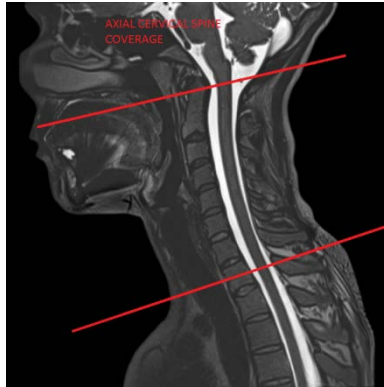

AXIAL Coverage