


Cervical Spine Trauma

Application Examples: trauma

Sequence	FoV	Thickness	Gap
SAG T1	24CM	3MM	0.5MM
SAG T2	24CM	3MM	0.5MM
SAG STIR	24CM	3MM	0.5MM
SAG 2D MERGE / MEDIC (SIEMENS)	24CM	3MM	0.5MM
SAG PD	24CM	3MM	0.5MM
AX 2D MERGE / MEDIC (SIEMENS)	18CM	3MM	0.3MM
AX T2	18CM	3MM	0.3MM



SAG Coverage

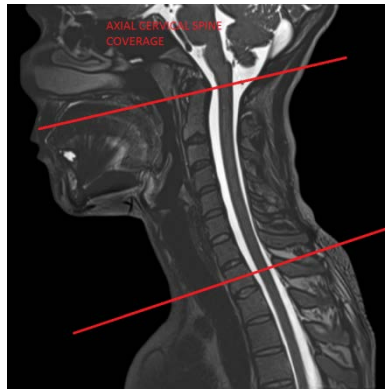

AXIAL Coverage