


## Thoracic Spine with and without GAD

Application Examples: history of cancer, infection, MS

Sequence	FoV	Thickness	Gap
SAG T1	35CM	3MM	0.5MM
SAG T2	35CM	3MM	0.5MM
SAG STIR	35CM	3MM	0.5MM
AX T2 UPPER <i>(MAKE SURE TO COVER ENTIRE TSPINE)</i>	20CM	4MM	0.5MM
AX T2 LOWER <i>(MAKE SURE TO COVER ENTIRE TSPINE)</i>	20CM	4MM	0.5MM
SAG DWI <i>(VALUES B-0,B-1000)</i>	35CM	3MM	1MM
INJECT CONTRAST (1mmol/Kg)			
SAG T1 FS WITH GAD	35CM	3MM	0.5MM
AX T1 FS UPPER WITH GAD <i>(MAKE SURE TO COVER ENTIRE TSPINE)</i>	20CM	4MM	0.5MM
AX T1 FS LOWER WITH GAD <i>(MAKE SURE TO COVER ENTIRE TSPINE)</i>	20CM	4MM	0.5MM



SAG Coverage

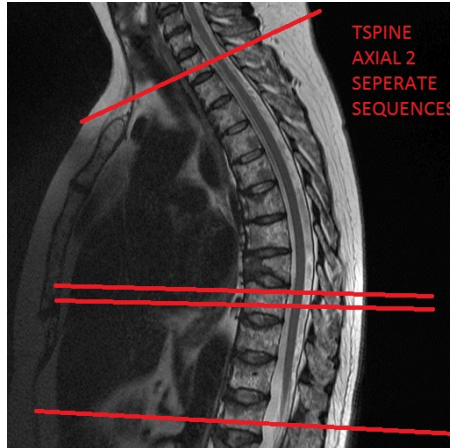

AXIAL Coverage