



Lumbosacral Plexus (Sciatic Neurogram) 3T ONLY

Application Examples: sciatic nerve pain

Sequence	FoV	Thickness	Gap
AX T1	40CM	3MM	1MM
AX T2 FS	40CM	3MM	1MM
COR STIR	40CM	3MM	1MM
COR T1	40CM	3MM	1MM
INJECT CONTRAST (1mmol/kg)			
AX T1 FS WITH GAD	40CM	3MM	1MM
COR T1 FS WITH GAD	40CM	3MM	1MM

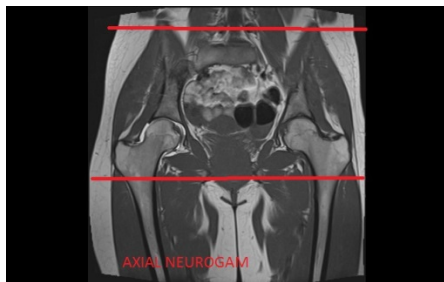

AXIAL Coverage

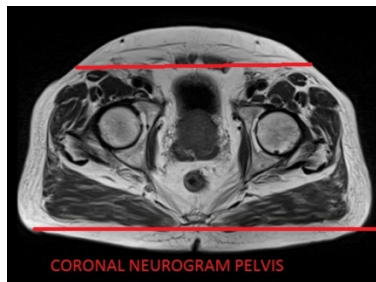

COR Coverage