



## MRV Head

Application Examples: thrombosis, headache, papilledema, hypertension

Sequence	FoV	Thickness	Gap
AX SWAN SWI	24CM	3MM	NO GAP
SAG 3D MRV TOF	24CM	1MM	NO GAP
AX T1 FS	24CM	4MM	1MM
AX BRAVO PRE	24CM	2MM	NO GAP
INJECT CONTRAST 1mmol/kg			
AX BRAVO WITH GAD ( <i>SUBTRACT PRE AND POST, SEND SUBTRACTED IMAGES ONLY</i> )	24CM	2MM	NO GAP

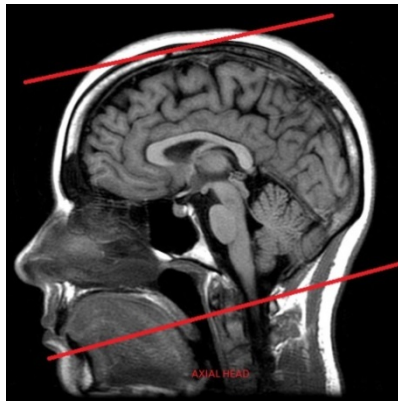

AXIAL Coverage

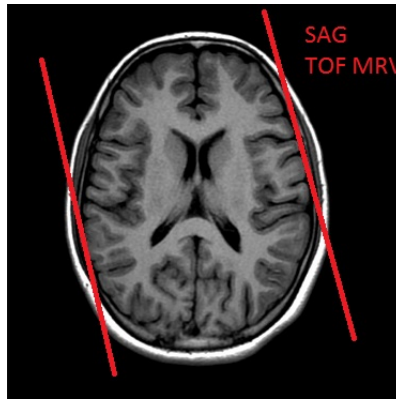

SAG Coverage