


MRA Head

Application Examples: family history of aneurysm, post coil, headaches

Sequence	FoV	Thickness	Gap
AX 3D TOF PRE **RECON MIPS, COR AND SAG RECON OF SOURCE IMAGES**	22CM	1.4MM	NO GAP
INJECT CONTRAST 1mmol/kg			
COR CE MRA FT (ONLY IF CONTRAST IS ORDERED) **RECON MIPS, AX AND SAG RECON OF SOURCE IMAGES**	26CM	1.8MM	NO GAP

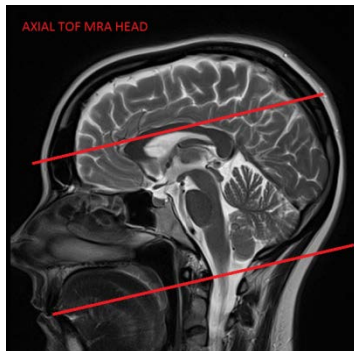

AXIAL Coverage



MIP