

Routine Ankle

Application Examples: Pain, Spain, PTT Dysfunction

Sequence	FoV	Thickness	Gap
Axial T1	14.0	3.0mm	0.8mm
Axial T2 FS	14.0	3.0mm	0.8mm
Coronal T1	14.0	3.0mm	0.8mm
Coronal T2 FS	14.0	3.0mm	0.8mm
Sagittal T1	14.0	3.0mm	0.8mm
Sagittal T2 FS	14.0	3.0mm	0.8mm

Anatomy Coverage: base of 5th Metatarsal through entire Calcaneus and Achilles tendon

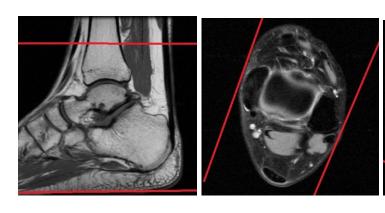

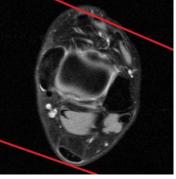

Axial Coverage

Sagittal Coverage

Coronal Coverage