


Routine Foot


Application Examples: Pain, fracture

Sequence	FoV	Thickness	Gap
Axial T1	14.0	3.0mm	0.8mm
Axial T2 FS	14.0	3.0mm	0.8mm
Coronal T1	14.0	3.0mm	0.8mm
Coronal T2 FS	14.0	3.0mm	0.8mm
Sagittal T1	14.0	3.0mm	0.8mm
Sagittal T2 FS	14.0	3.0mm	0.8mm

Anatomy Coverage: Navicular though toes

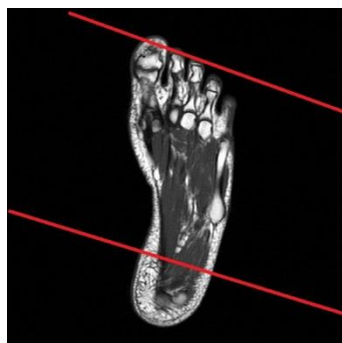

Axial Coverage

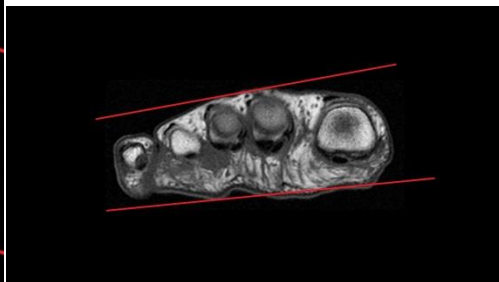

Coronal Coverage



Sagittal Coverage