


Routine Knee


Application Examples: Pain, meniscal injury, ACL injury, internal derangement, Osteoarthritis

If patient has had prior surgery and has metal in knee please run Stir sequences instead of T2 FS

Sequence	FoV	Thickness	Gap
Axial T2 FS/Stir	15.0	3.0mm	1mm
Coronal PD	15.0	3.0mm	0.5mm
Coronal T2 FS/Stir	15.0	3.0mm	0.5mm
Sagittal PD	15.0	3.0mm	0.5mm
Sagittal T2 FS/Stir	15.0	3.0mm	0.5mm

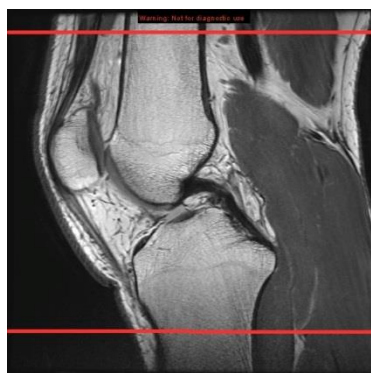

Axial Coverage

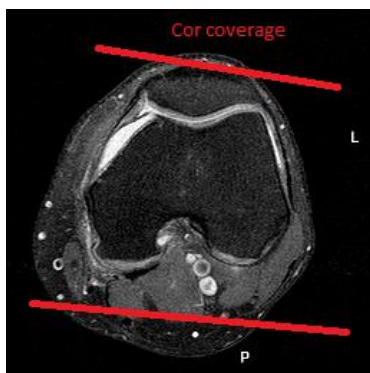

Coronal Coverage



Sagittal Coverage