


Routine Hip Arthrogram


Application Examples: Laberal Tear

All images should be straight to body


Sequence	FoV	Thickness	Gap
Axial T2 FS entire Pelvis	35.0	4.0mm	1.0mm
Axial Fiesta (<i>Recon 1x1 in both Sagittal and Coronal plane</i>)	25.0	1.0mm	No Gap
Axial T1 Non FS	18.0	3.0mm	0.5mm
Coronal T1 FS	18.0	3.0mm	0.5mm
Coronal T2 FS	18.0	3.0mm	0.5mm
Sagittal T1 FS	18.0	3.0mm	0.5mm
Axial T1 Obl FS	18.0	3.0mm	0.5mm

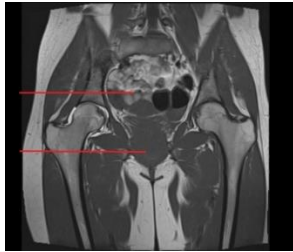


Axial Coverage

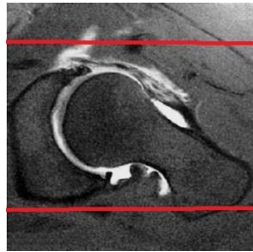

Coronal Coverage

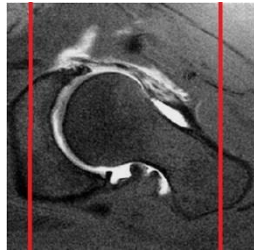

Sagittal Coverage

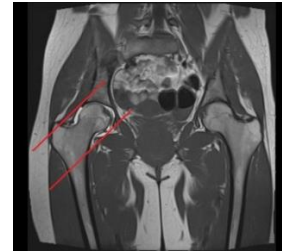

Axial Obl Coverage