



Routine Wrist

Application Examples: Wrist pain, scaphoid fx, ganglion cyst


Sequence	FoV	Thickness	Gap
Axial T1	10.0	2.5mm	0.5mm
Axial T2 FS	10.0	2.5mm	0.5mm
Coronal PD	10.0	2.5mm	0.5mm
Coronal T2 FS	10.0	2.5mm	0.5mm
Coronal 3D Merge	10.0	2.5mm	0.5mm
Sagittal T2 FS	10.0	2.5mm	0.5mm

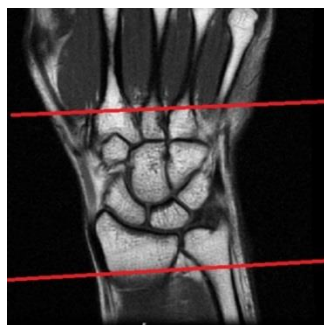

Axial Coverage

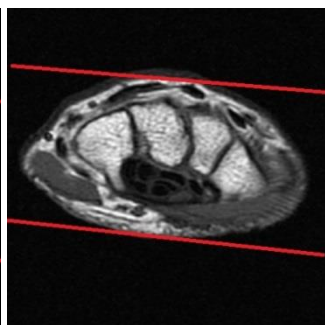

Coronal Coverage

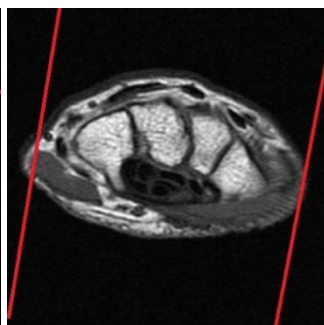

Sagittal Coverage