



Routine Wrist Arthrogram

Application Examples: TFCC tear, loose body


Sequence	FoV	Thickness	Gap
Axial T1 FS	10.0	2.5mm	0.5mm
Axial T2 FS	10.0	2.5mm	0.5mm
Coronal T1 FS	10.0	2.5mm	0.5mm
Coronal T1	10.0	2.5mm	0.5mm
Coronal T2 FS	10.0	2.5mm	0.5mm
Sagittal T1	10.0	2.5mm	0.5mm
Sagittal T2	10.0	2.5mm	0.5mm

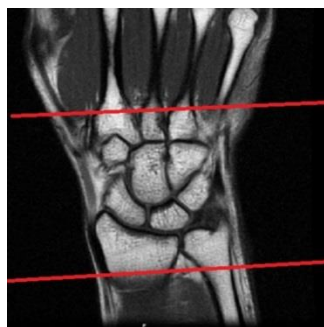

Axial Coverage

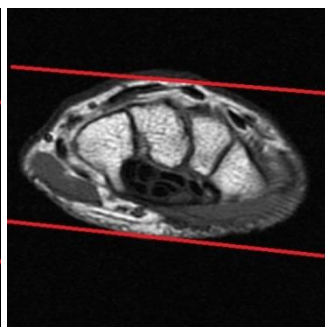

Coronal Coverage

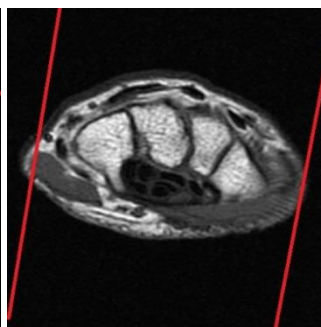

Sagittal Coverage